


Concerned about the threat of antibiotic resistance?



Do you expect a precise identification
of your carbapenemase?

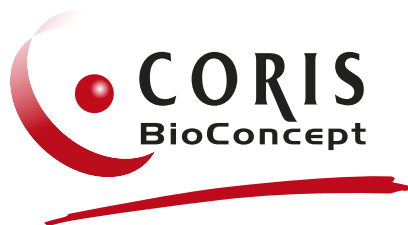

Innovative solutions for effective diagnostics

Antibio-resistance is rising up !

Carbapenemase-producing Organisms (CPO) represent a major global health concern due to their wide spectrum of resistance to most classes of antimicrobial agents, thus leaving very few remaining options for the therapeutic management of infected patients. Due to their extremely fast spreading, CPO needs to be rapidly identified in order to isolate harbouring patients and to adopt the most accurate antibiotic treatment.

Are you aware about the consequences ?




CPO confirmation

Isolation of carriers


Surveillance of close contacts

Healthcare unit on stand by

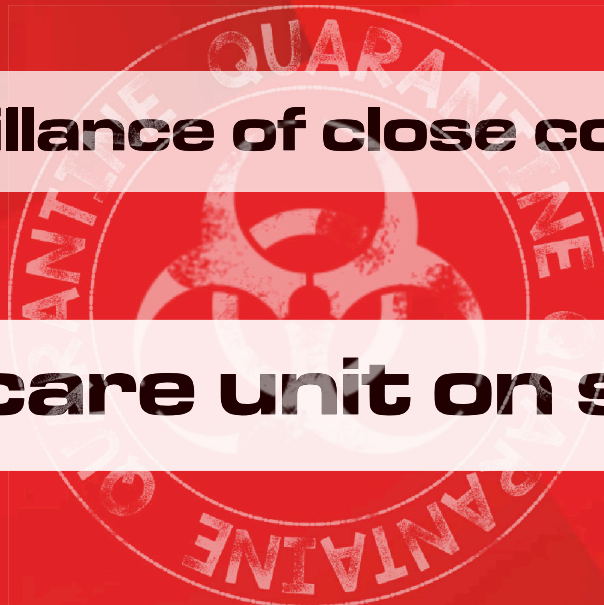


Tackle the three main carbapenemases in Europe !

OXA-48

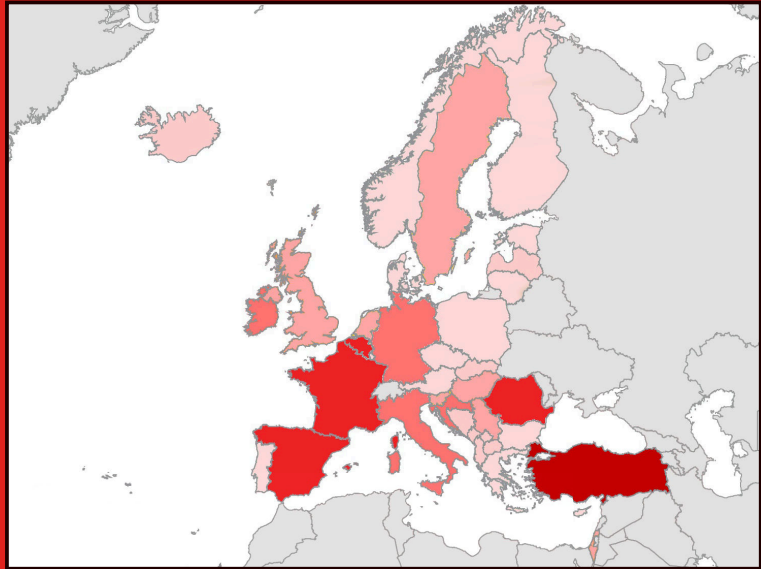


© EUOROUNDUP : Carbapenemase-producing Enterobacteriaceae in Europe: assessment by national experts from 38 countries, May 2015 - B. Albiger, C. Glasner, M.J. Struelens, H. Grundmann, D.L. Monnet, the European Survey of Carbapenemase-Producing Enterobacteriaceae (EuSCAPE) working group

KPC

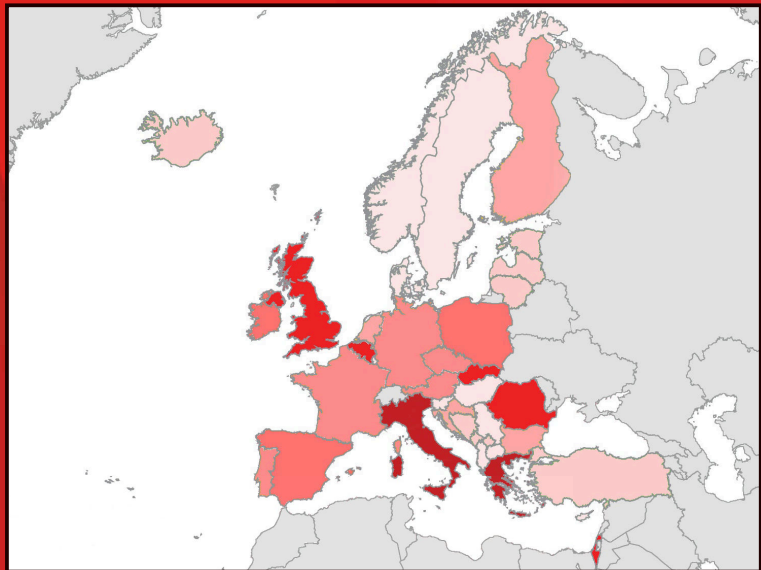

© EUOROUNDUP : Carbapenemase-producing Enterobacteriaceae in Europe: assessment by national experts from 38 countries, May 2015 - B. Albiger, C. Glasner, M.J. Struelens, H. Grundmann, D.L. Monnet, the European Survey of Carbapenemase-Producing Enterobacteriaceae (EuSCAPE) working group

NDM

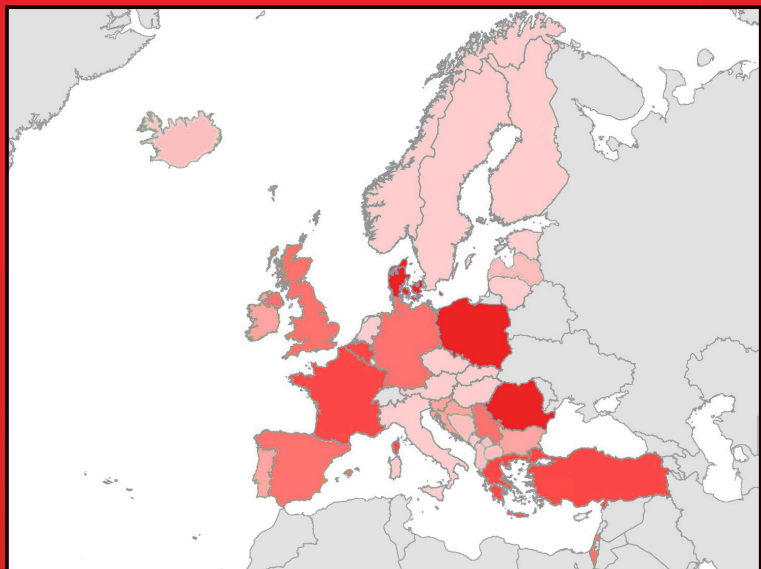

© EUOROUNDUP : Carbapenemase-producing Enterobacteriaceae in Europe: assessment by national experts from 38 countries, May 2015 - B. Albiger, C. Glasner, M.J. Struelens, H. Grundmann, D.L. Monnet, the European Survey of Carbapenemase-Producing Enterobacteriaceae (EuSCAPE) working group

A unique solution !

**Coris BioConcept offers a unique test to identify
3 Carbapenemases in only 1 sample !**



RESIST-3 O.K.N.

- 3 parameters in 1 single test
- No equipment required
- Easy and ready to use
- Shortest Time-To-Result
- High performances

RESIST-3 O.K.N.

OXA-48

KPC


NDM

Technical details

Carbapenems are β -lactam antibiotics mainly used in hospital for the treatment of nosocomial multidrug resistant bacteria.

Carbapenemases are bacterial enzymes that hydrolyse carbapenem molecules, providing antibiotic-resistance to Carbapenemase-producing Organisms (CPO). Resistances conferred by the 5 main carbapenemases (OXA-48, KPC, NDM, IMP and VIM) continue to increase. Any diagnostic tool allowing to help their laboratory identification can improve patient's safety.

In its commitment to offer adapted diagnostic solutions, Coris BioConcept recently developed a new family of immunochromatography tests under the brand name "RESIST". This innovative range of rapid antigen tests is unique worldwide. Due to their high performances, RESIST tests have already been adopted by laboratories for a fast and reliable identification of antibiotic-resistance bacteria directly from culture colonies.

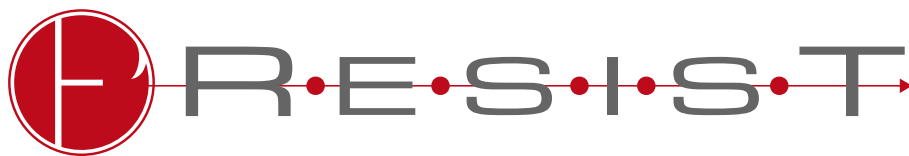

Performances

OXA-48

Molecular method (n=173)		
Specificity	100%	(97.2 to 100%)
Sensitivity	100%	(95.7 to 100%)
PPV	100%	(95.7 to 100%)
NPV	100%	(97.2 to 100%)
Reliability	100%	


KPC

Molecular method (n=173)		
Specificity	100%	(98.2 to 100%)
Sensitivity	100%	(68.4 to 100%)
PPV	100%	(68.4 to 100%)
NPV	100%	(98.2 to 100%)
Reliability	100%	

NDM

Molecular method (n=173)		
Specificity	100%	(98.1 to 100%)
Sensitivity	100%	(84.7 to 100%)
PPV	100%	(84.7 to 100%)
NPV	100%	(98.1 to 100%)
Reliability	100%	

Principle of use

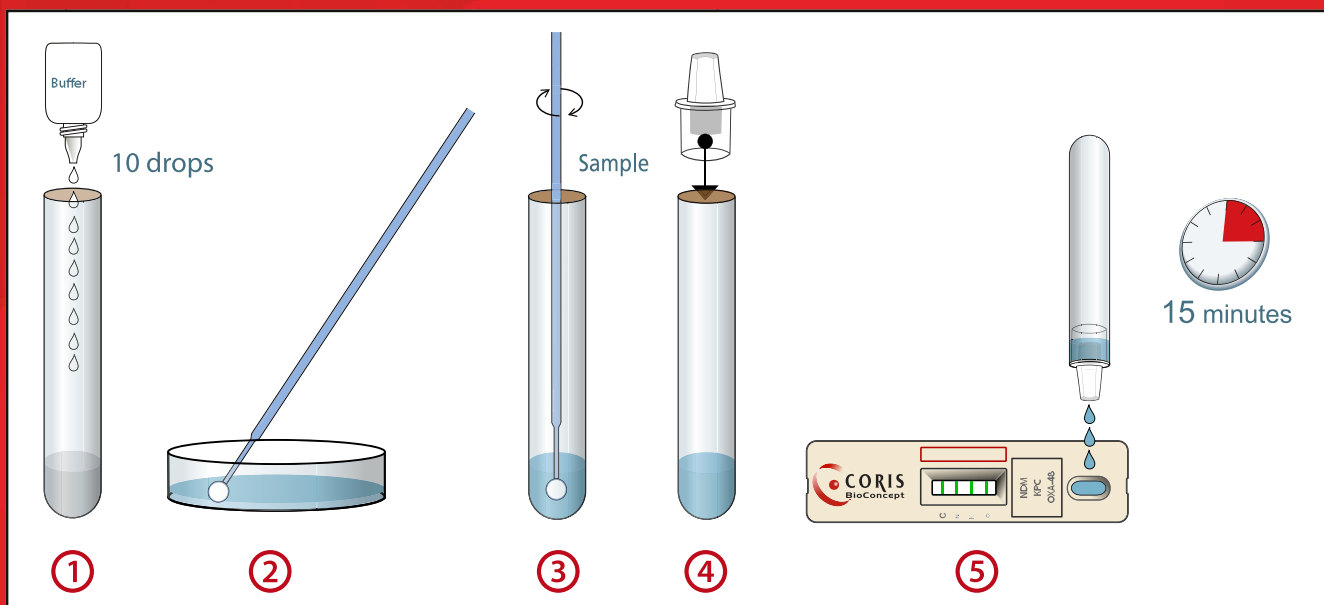

Ordering informations

Product name	Target	Description	Storage	Code
RESIST-3 O.K.N. K-SeT	OXA-48-like KPC-like NDM-like	<ul style="list-style-type: none">• 20 individually sealed pouches• 15 mL dilution buffer in a bottle• 20 collection tubes with droppers• Instruction For Use	4°C to 30°C	K-15R5

Also available:

OXA-48 K-SeT	OXA-48-like	<ul style="list-style-type: none">• 20 individually sealed pouches• 15 mL dilution buffer in a bottle• 20 collection tubes with droppers• Instruction For Use	4°C to 30°C	K-15R1
KPC K-SeT	KPC-like	<ul style="list-style-type: none">• 20 individually sealed pouches• 15 mL dilution buffer in a bottle• 20 collection tubes with droppers• Instruction For Use	4°C to 30°C	K-15R2

- Price: please refer to your local distributor
- Stability: up to 12 months
- Bibliographic and Technical support: client.care@corisbio.com

Distributed by:



www.corisbio.com

Coris BioConcept Headquarters

Science Park CREALYS - Rue Jean Sonet 4A
5032 Gembloux - Belgium
Tel: +32 (0)81 719 917
Fax: +32 (0)81 719 919
E-mail: info@corisbio.com - sales@corisbio.com
www.corisbio.com

Coris BioConcept R&D Department

Science Park CREALYS - Rue Jean Sonet 29, Bte 28
5032 Gembloux - Belgium
Tel: +32 (0)81 559 320
Fax: +32 (0)81 719 919
E-mail: info@corisbio.com - sales@corisbio.com
www.corisbio.com